

anti-human CD4 FITC-conjugated

FITC - conjugated monoclonal antibody to CD4 (Human)

Cat-No: **21850043**

500 µl

Clone: OKT-4

Specificity: The OKT4 monoclonal antibody reacts with human CD4, a 59 kDa cell surface glycoprotein expressed by the majority of thymocytes, a subpopulation of mature T cells (T-helper cells) and in low levels on monocytes. CD4 is a receptor for the human immunodeficiency virus (HIV). The OKT4 antibody recognizes a different epitope than the RPA-T4 monoclonal antibody, and these antibodies do not cross-block binding to each other's respective epitopes.

Isotype subclass: Mouse IgG2b

Physical state: Liquid

Buffer/Additives/Preservative: PBS containing 1 % BSA and 0.09 % sodium azide (pH 7.2)

Form: The purified antibody is conjugated with Fluorescein isothiocyanate (FITC) under optimum conditions. The reagent is adjusted for direct use. No reconstitution is necessary.

Expiration date: The reagent is stable until the expiry date stated on the vial label.

Storage conditions: Store at 4 °C. Avoid prolonged exposure to light.


Application: Flow cytometry

References:

- 1) Bour S, Boulerice F, Wainberg MA. Inhibition of gp160 and CD4 maturation in U937 cells after both defective and productive infections by human immunodeficiency virus type 1. J Virol. 1991 Dec;65(12):6387-96. (OKT4, IP, PubMed)
- 2) Reinherz EL, Kung PC, Goldstein G, Schlossman SF. 1979. Separation of functional subsets of human T cells by a monoclonal antibody. Proc Natl Acad Sci. 76(8): 4061-5.

Warning: Sodium azide is harmful if swallowed (R22). Keep out of reach of children (S2). Keep away from food, drink and animal feeding stuff (S13). Wear suitable protective clothing (S36). If swallowed, seek medical advice immediately and show this container or label (S46). Contact with acids liberates very toxic gas (R32). Azide compounds should be flushed with large volumes of water during disposal to avoid deposits in lead or copper plumbing where explosive conditions can develop.

This material is offered for **research use only**. Not for use in human. For in vitro use only. ImmunoTools will not be held responsible for patent infringement or other violations that may occur with the use of our products.

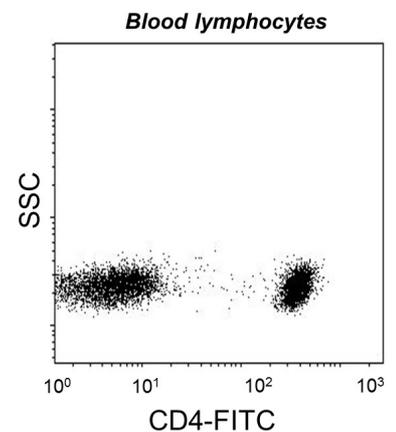

Attention! Cells from one healthy individual are shown. Cell Populations and staining intensity may vary interindividually.

ImmunoTools Excellent Quality - Advantageously priced

Gladiolenweg 2; 26169 Friesoythe; Germany

phone:+49-(0)4491-400997, fax:+49-(0)4491-400998, info@immunotools.com

www.immunotools.com